

# The Eliminator 200 Series

The Eliminator 200 series provides the same aerodynamic design features as the EH-100 airfoil hood series, with entrance from two sides. This makes it ideal for use in classroom or prep room settings or whenever demonstration of experiments is desired. A 45-degree angle around the fascia and a flush-mount radiused airfoil across the bottom provide minimized turbulence and increased performance.


## ELIMINATOR 200 - DOUBLE-SIDED 100 SERIES

### Double-Sided 100 Series - Dual Access



**U.L. 1805  
Classified**

- Powder coat finish
- Full frame construction
- U.L. listed poly resin liner
- Flush-mount airfoil
- Knock down capability
- LED light with bulbs
- Full-length finger lift
- Louvered front for by-pass
- 10" round S/S duct collar

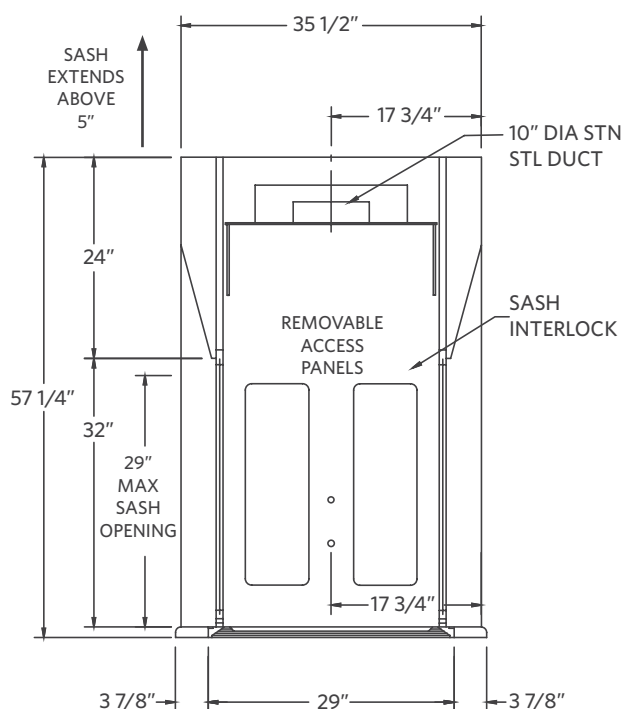






| Part Number | Superstructure |     |         | Worktop |        |     | Hood Opening |     | CFM @ 29"* | CFM @ 18"* | SP  | Blower Model # | Shipping Weight |
|-------------|----------------|-----|---------|---------|--------|-----|--------------|-----|------------|------------|-----|----------------|-----------------|
|             | W              | H   | D       | W       | H      | D   | W            | H   |            |            |     |                |                 |
| EH-211-36   | 36"            | 56" | 35-1/2" | 36"     | 1-1/4" | 29" | 26"          | 29" | 567        | 354        | .25 | BL-100-36      | 400             |
| EH-211-48   | 48"            | 56" | 35-1/2" | 48"     | 1-1/4" | 29" | 38"          | 29" | 833        | 521        | .33 | BL-100-48      | 450             |
| EH-211-60   | 60"            | 56" | 35-1/2" | 60"     | 1-1/4" | 29" | 50"          | 29" | 1100       | 688        | .34 | BL-100-60      | 500             |
| EH-211-72   | 72"            | 56" | 35-1/2" | 72"     | 1-1/4" | 29" | 62"          | 29" | 1367       | 854        | .34 | BL-100-72      | 550             |
| EH-211-96   | 96"            | 56" | 35-1/2" | 96"     | 1-1/4" | 29" | 86"          | 29" | 1900       | 1188       | .40 | BL-100-96      | 650             |

# The Eliminator 200 Series

The EH-200 features standard 5" (127 mm) posts to allow for electrical and remote control plumbing services on both sides, and these can be configured to your exact specifications.

Five standard sizes and an array of options make this one of the most flexible fume hoods available anywhere.

## ELIMINATOR 200 - DOUBLE-SIDED 100 SERIES

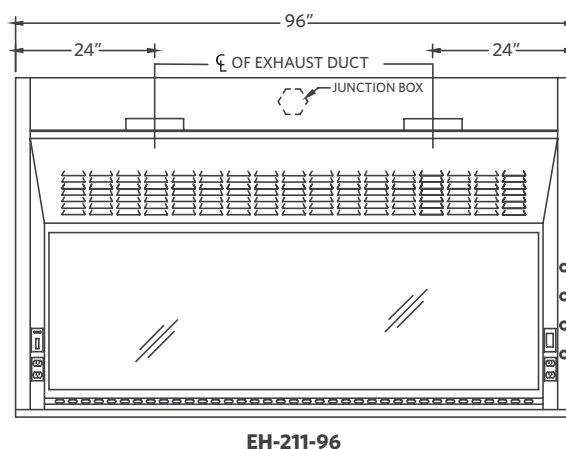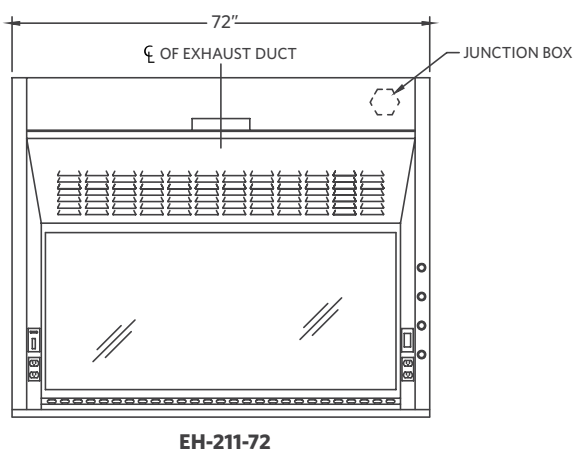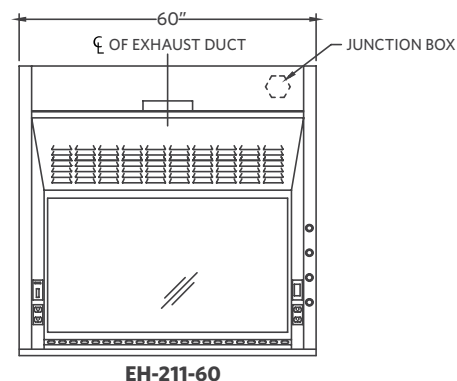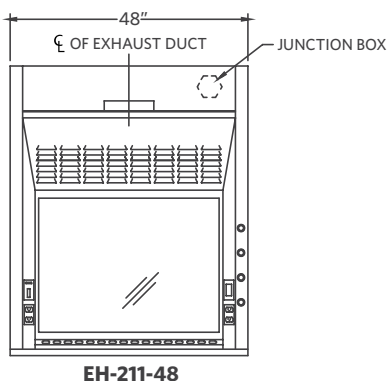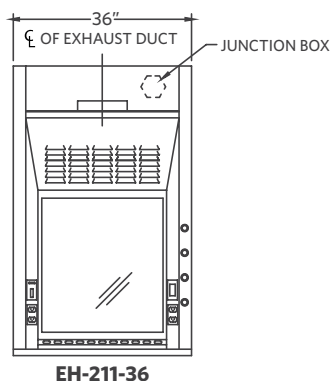

## Accessories

| ITEM                   | AMS PT# | PAGE # | ITEM                     | AMS PT#     | PAGE # |
|------------------------|---------|--------|--------------------------|-------------|--------|
| Ceiling Enclosures     | CE      | 33     | Sash Stops               | SASHSTOP    | 38     |
| Plumbing Accessories   |         | 34     | Chain and Sprocket Drive | CHNSPK      | CALL   |
| Electrical Accessories |         | 35     | Stainless Steel Sash     | SS SASH     | 38     |
| Epoxy Tops             | ET      | 39     | Stainless Steel Airfoil  | SS AF       |        |
| Alarms                 | ALARM   | 32     | Push Button Sash         | PUSH BUTTON | 38     |
| Base Cabinets          | BC      | 36     | Combination Sash         | H&V         | 38     |
| Blowers                | BL      | CALL   | Lattice Assembly         | LATTICE     | CALL   |
| Pre-Plumbing           | PREP    |        | Auto Sash Return         | AUTO RETURN | 38     |
| Pre-Wiring             | PREW    |        |                          |             |        |